

## **IC15 HOMECARE MATTRESS**

The IC15 offers support and comfort. Firm feel makes it easier for turnina.

Pair with an Icare bed or suitable on any base platform.

## **FIRM FEEL**





#### Pressure Redistribution

The IC15 mattress features the ActiveX 84 material is the main comfort layer in this mattress. ActiveX is a high grade of elastic type foam which is responsive to temperature around inflamed areas. ActiveX will become softer as it becomes warmer. So around a pressure point it will always soften and allow the blood to circulate freely.



### Solid Support Foam Base

XP Support foam is a highdensity foam that is used in the bottom layer of the IC15. XP Support foam is an open cell material that breaths to allow air circulation but is high density to prevent bottoming out. This material will easily contour to any profiling bed.



### **Temperature Regulating Cover**

The IC15 cover is made with the new Bluetec cooling fabric. This is all way stretch fabric to conform with body and mattress contour. It is under woven with a cooling yarn to regulate body temperature.

KING SINGLE

Code: IC15KS



### Ventilated Base Fabric Vant fabric is a three-

dimensional polypropylene material that has high air permeability and is antibacterial. We use this on the side wall of the IC15 mattress covers to greatly increase air movement.



#### Medical Covers Available

Zerotec fabric is a medical barrier fabric that's has high elasticity in every direction which reduces shear on pressure redistribution surfaces. It is waterproof and wipeable but has microporous breathable properties. This cover has welded seams and waterfall flaps for fluid control.

## CODES & SIZES:



## Mattress Size: 2030 x 900 QUEEN

LONG SINGLE

Code: IC15LS

Code: IC15Q Mattress Size: 2030 x 1520



## **SPLIT QUEEN** Code: IC15SQ

Mattress Size: 2030 x 1070

Mattress Size: 2030 x 760 (Per side)



### **SPECIFICATIONS**

- No turn design reducing manual handling »
- Anticrease XP base to allow easy contour for any profiling bed »
- » Standard cover MVTR: NA. Medical Cover MVTR (g/24h·m2): 1635.0 Average
- Medical Cover BS7175: Crib 5 »
- Waterlow ≤20 »
- To launder the cover in a washing machine, set the main wash to a minimum of » 65°C for no less than 10 minutes. The cover can be tumble dried on low heat and the temperature must not exceed 60°C. The foam can be autoclaved at 134°C.











## **HIGH RISK**

This rating is based on our average peak pressure of <24mmHg. For some high risk clients Icare Zerotec covers are recommended.

# **IC20 HOMECARE MATTRESS**



Middle of the range, the IC20 proves to be popular for most. The medium feel provides a soft feel but excellent support.

Medical cover's available if required.

## **MEDIUM FEEL**





### Pressure Redistribution

The IC20 mattress features the ActiveX 84 material is the main comfort layer in this mattress. ActiveX is a high grade of elastic type foam which is responsive to temperature around inflamed areas. ActiveX will become softer as it becomes warmer. So around a pressure point it will always soften and allow the blood to circulate freely.



### Solid Support Foam Base

XP Support foam is a highdensity foam that is used in the bottom layer of the IC20. XP Support foam is an open cell material that breaths to allow air circulation but is high density to prevent bottoming out. This material will easily contour to any profiling bed.



### **Temperature Regulating Cover**

The IC20 cover is made with the new Bluetec cooling fabric. This is all way stretch fabric to conform with body and mattress contour. It is under woven with a cooling yarn to regulate body temperature.

LONG SINGLE

Code: IC2LS

LONG DOUBLE

Code: IC2LD

Mattress Size: 2030 x 900

Mattress Size: 2030 x 1350



## Ventilated Base Fabric

Vant fabric is a threedimensional polypropylene material that has high air permeability and is antibacterial. We use this on the side wall of the IC20 mattress covers to greatly increase air movement.



#### Medical Covers Available

Zerotec fabric is a medical barrier fabric that's has high elasticity in every direction which reduces shear on pressure redistribution surfaces. It is waterproof and wipeable but has microporous breathable properties. This cover has welded seams and waterfall flaps for fluid control.

### CODES & SIZES:



## SINGLE

Code: IC2S Mattress Size: 1900 x 900



## QUEEN

Code: IC2Q Mattress Size: 2030 x 1520



## SPLIT QUEEN Code: IC2SQ


Mattress Size: 2030 x 760 (Per side)

## KING SINGLE Code: IC2KS Mattress Size: 2030 x 1070

SPECIFICATIONS

- » No turn design reducing manual handling
- Anticrease XP base to allow easy contour for any profiling bed »
- Standard cover MVTR: NA. Medical Cover MVTR (g/24h·m2): 1635.0 Average »
- Medical Cover BS7175: Crib 5 »
- Waterlow ≤20 »
- To launder the cover in a washing machine, set the main wash to a minimum of 65°C for no >> less than 10 minutes. The cover can be tumble dried on low heat and the temperature must not exceed 60°C. The foam can be autoclaved at 134°C.



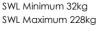







**HIGH RISK** 

This rating is based on our average peak pressure of <19mmHg. For some high risk clients Icare Zerotec covers are recommended.



## **IC25 HOMECARE MATTRESS**

For clients who like a soft mattress, or who spend a long time in bed. The feel of the IC25 Mattress is best suited to aged or petite patients.



**SOFT FEEL** 



#### Pressure Redistribution

The IC25 mattress features the ActiveX 84 material is the main comfort layer in this mattress. ActiveX is a high grade of elastic type foam which is responsive to temperature around inflamed areas. ActiveX will become softer as it becomes warmer. So around a pressure point it will always soften and allow the blood to circulate freely.



### Solid Support Foam Base

XP Support foam is a highdensity foam that is used in the bottom layer of the IC25. XP Support foam is an open cell material that breaths to allow air circulation but is high density to prevent bottoming out. This material will easily contour to any profiling bed.



### **Temperature Regulating Cover**

The IC25 cover is made with the new Bluetec cooling fabric. This is all way stretch fabric to conform with body and mattress contour. It is under woven with a cooling yarn to regulate body temperature.



### Ventilated Base Fabric Vant fabric is a three-

dimensional polypropylene material that has high air permeability and is antibacterial. We use this on the side wall of the IC25 mattress covers to greatly increase air movement.



### Medical Covers Available

Zerotec fabric is a medical barrier fabric that's has high elasticity in every direction which reduces shear on pressure redistribution surfaces. It is waterproof and wipeable but has microporous breathable properties. This cover has welded seams and waterfall flaps for fluid control.

## CODES & SIZES:



SINGLE Code: ICXLS Mattress Size: 1900 x 900



## QUEEN

Code: ICXLQ Mattress Size: 2030 x 1520





## LONG SINGLE Code: ICXLLS Mattress Size: 2030 x 900





# Mattress Size: 2030 x 1350

LONG DOUBLE

Code: ICXLLD

SPLIT QUEEN Code: ICXLSQ Mattress Size: 2030 x 760 (Per side)

### SPECIFICATIONS

- No turn design reducing manual handling »
- Anticrease XP base to allow easy contour for any profiling bed »
- Standard cover MVTR: NA. Medical Cover MVTR (g/24h·m2): 1635.0 Average »
- Medical Cover BS7175: Crib 5 »
- » Waterlow ≤20
- 23 To launder the cover in a washing machine, set the main wash to a minimum of 65°C for no less than 10 minutes. The cover can be tumble dried on low heat and the temperature must not exceed 60°C. The foam can be autoclaved at 134°C.



SWL Minimum 32kg SWL Maximum 212kg







## **HIGH RISK**

This rating is based on our average peak pressure of <22mmHg. For some high risk clients Icare Zerotec covers are recommended.

Icare risk category ratings are based on average peak pressure, Therapist recommendation, case studies, industry comparison and the Waterlow scale. For some high risk clients Icare Zerotec covers are recommended. The best way to assess a client on a pressure reducing surface is to monitor pressure points with visual checks over a period of 5 days.


ICARE MEDICAL GROUP - HOMECARE BEDS & MATTRESSES CATALOGUE

## **M2 MEDICAL MATTRESS**



The M2 mattress features firm edge support and high density foams. This is a high risk mattress, and suitable for the use on all Icare beds.

### **MEDIUM FEEL**

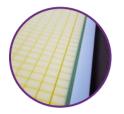

#### Pressure Redistribution

XP Light foam is very responsive and microporous material allowing a high level of airflow. This is used in the top layer of the M2 mattress and is module cut to allow each square to move to body contour. This module cutting also allows for horizontal anti shear movement.



## Base & Rail Support Foam

XP Support foam is a highdensity foam that is used in the bottom layer and side rails of the M2. XP Support foam is an open cell material that breaths to allow air circulation but is high density to prevent bottoming out or side roll out. This material will easily contour to any profiling bed.



### Zerotec Cover

Zerotec fabric is a medical barrier fabric that has high elasticity in every direction which reduces shear on pressure redistribution surfaces. It is waterproof and wipeable but has microporous breathable properties. This is used for the cover of the M2 mattress. This cover has welded seams and waterfall flaps for fluid control.

## CODES & SIZES:



LONG SINGLE Code: ICM2LS Mattress Size: 2030 x 900

LONG DOUBLE Code: ICM2LD Mattress Size: 2030 x 1350



**KING SINGLE** Code: ICM2KS Mattress Size: 2030 x 1070 QUEEN

Code: ICM2Q Mattress Size: 2030 x 1520

### SPECIFICATIONS

- No turn design reducing manual handling »
- Anticrease XP Support base to allow easy contour for any profiling bed »
- Standard cover MVTR: NA. Medical Cover MVTR (g/24h·m2): 1635.0 Average »
- Medical Cover BS7175: Crib 5 >>
- Waterlow ≤20 »
- )) To launder the cover in a washing machine, set the main wash to a minimum of 65°C for no less than 10 minutes. The cover can be tumble dried on low heat and the temperature must not exceed 60°C. The foam can be autoclaved at 134°C.











## **HIGH RISK**


This rating is based on our average peak pressure of <20mmHg. For more information refer to the Icare pressure care manual.



## **M4 MEDICAL MATTRESS**

The M4 mattress is very high risk profiling foam mattress built with a high density base. The mattress is designed with a section of v-cut XP Support Foam which creates additional air movement for enhanced pressure relieving properties.

## **MEDIUM FEEL**

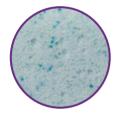

### Pressure Redistribution

ActiveGEL material is a high grade of elastic type foam infused with cooling gel partials. This material is responsive to temperature and will become softer as it becomes warmer. So around a pressure point it will always soften and allow the blood to circulate freely. The Gel infusion helps to regulate overall temperature.



Castellated ActiveGEL The ActiveGEL is module cut to allow each square to move to body contour. This module cutting also allows for horizontal anti shear

movement.



### Base & Rail Support Foam

XP Support base layer features deep cut springs allowing much more airflow and extension of lateral contouring. XP Support foam is an open cell material that

breaths to allow air circulation but is high density to prevent bottoming out or side roll out. This material will easily contour to any profiling bed.

**KING SINGLE** 



Zerotec Cover

Zerotek fabric is a medical barrier fabric that has high elasticity in every direction which reduces shear on pressure redistribution surfaces. It is waterproof and wipeable but has microporous breathable properties. This is used for the cover of the M2 mattress. This cover has welded seams and waterfall flaps for fluid control.

### CODES & SIZES:



### LONG SINGLE Code: ICM4LS Mattress Size: 2030 x 900

## LONG DOUBLE

Code: ICM4LD Mattress Size: 2030 x 1350



Code: ICM4KS Mattress Size: 2030 x 1070 **QUEEN** 

Code: ICM4Q Mattress Size: 2030 x 1520

### SPECIFICATIONS

- » No turn design reducing manual handling
- » Anticrease XP Support base to allow easy contour for any profiling bed
- » Standard cover MVTR: NA. Medical Cover MVTR (g/24h·m2): 1635.0 Average
- » Medical Cover BS7177: Crib 5
- » Waterlow ≤20
- » To launder the cover in a washing machine, set the main wash to a minimum of 65°C for no less than 10 minutes. The cover can be tumble dried on low heat and the temperature must not exceed 60°C. The foam can be autoclaved at 134°C.











## HIGH RISK

This rating is based on our average peak pressure of <20mmHg. For more information refer to the Icare pressure care manual.

## LATEX MATTRESS



The Latex mattresses are a low risk pressure mattress but provide a unique bounce back support which makes it easy for transfers. Suitable on any base platform.

## **FIRM FEEL**



Natural Latex

Latex is naturally antibacterial. It is a natural rubber product that provides excellent support and bounce back.



Solid Support Foam Base The Latex core has pin hole ventilated giving the more air movement helping to regulate temperature



Quilted Cover Super soft quilted cover giving the mattress a home look and feel

**KING SINGLE** 

## CODES & SIZES:



LONG SINGLE Code: ICEFLS Mattress Size: 2030 x 900





Code: ICEFKS Mattress Size: 2030 x 1070



### SPECIFICATIONS

- No turn design reducing manual handling »
- Waterlow ≤14 »
- Latex material will easier contour to any profiling platform »
- Polycotton cotton tack and jump quilted cover »











## LOW RISK

This rating is based on our average peak pressure of <31mmHg. For some high risk clients Icare Zerotec covers are recommended.



## **MATTRESS OVERLAYS**

The range of Icare pressure care overlays are designed to provide added pressure relief and comfort on an existing mattress. These overlays are heat and pressure sensitive causing them to soften around pressure points and inflamed areas. These overlays come in all the mattress sizes up to queen.



## **MEDIUM FEEL**



### **Pressure Redistribution**

The overlay features the ActiveX material is the main comfort layer in this mattress. ActiveX is a high grade of elastic type foam which is responsive to temperature around inflamed areas. ActiveX will become softer as it becomes warmer. So around a pressure point it will always soften and allow the blood to circulate freely.



All-way Contour Cover The Icare overlay is made with an Elastex material cover. This fabric has a high elasticity in every direction enabling it to conform with body and

sleeping surface contour.



Ventilated Side Fabric Vant fabric is a threedimensional polypropylene material that has high air permeability and is antibacterial. We use this on the side wall of the overlay covers to greatly increase air movement.



Any Sleeping Surface

The Icare overlay can be used on any type of sleeping surface. The pressure redistribution and feel will vary depending on the under surface.

### CODES & SIZES:

SINGLE



Code: ICTS Mattress Size: 1900 x 900



DOUBLE Code: ICTD Mattress Size: 1900 x 1350



LONG SINGLE Code: ICTLS Mattress Size: 2030 x 900

Mattress Size: 2030 x 1520

QUEEN

Code: ICTQ



KING SINGLE Code: ICTKS Mattress Size: 2030 x 1070

### SPECIFICATIONS

- » Zip off washable cover
- » 50mm thickness
- » SWL range dependant on under surface foundation.
- » To launder the cover in a washing machine, set the main wash to a minimum of 65°C for no less than 10 minutes. The cover can be tumble dried on low heat and the temperature must not exceed 60°C. The foam can be autoclaved at 134°C.



SWL Minimum 32kg SWL Maximum 228kg



10 Years Warranty



HIGH RISK

This rating is based on our average peak pressure of <18mmHg. The pressure redistribution and feel will vary depending on the under surface. Some high risk clients may require a medical cover which is not available for this product.

Icare risk category ratings are based on average peak pressure, Therapist recommendation, case studies, industry comparison and the Waterlow scale. The best way to assess a client on a pressure reducing surface is to monitor pressure points with visual checks over a period of 5 days.